

Run Info

Experiment Name	ncov_ucdh_env1_run3
Sample ID	ncov_ucdh_env1_run3
Run ID	5aa17f3a-ca27-42a9-989b-61e1db6f36d6
Flow Cell Id	FAN19575
Start Time	September 21, 14:46
Run Length	23h 20m

Run Summary

Reads Generated	1.58 M
Bases Generated	306.08 Mb
Estimated Bases	523.44 Mb
Percentage Basecalled	63%


Run Parameters

Flow Cell Type	FLO-MIN106
Kit	SQK-LSK109
Basecalling	on
Specified Run Length	72 hours
Initial Bias Voltage	-180 mV
FAST5 Output	Enabled
FAST5 Output Options	zlib_compress,fastq,raw
FAST5 Reads per File	4000
FASTQ Output	Enabled
FASTQ Reads per File	4000
Active Channel Selection	Enabled
Mux Scan Period	1 hour 30 minutes
Reserved Pores	0 %
Basecall Model	Fast basecalling


Versions

MinKNOW Core	3.6.5
Bream	4.3.16
Guppy	3.2.10

Cumulative Output Reads

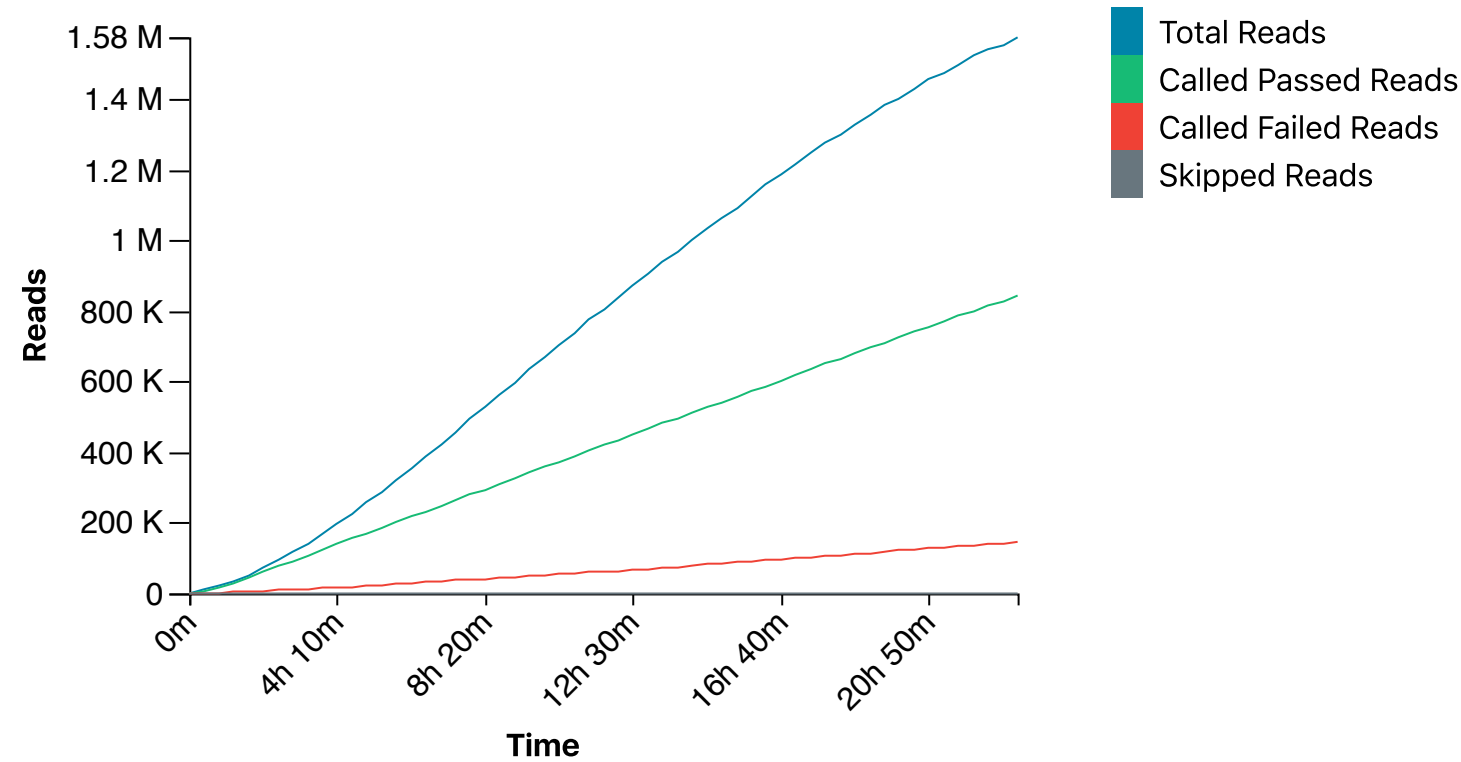


Cumulative Output Bases

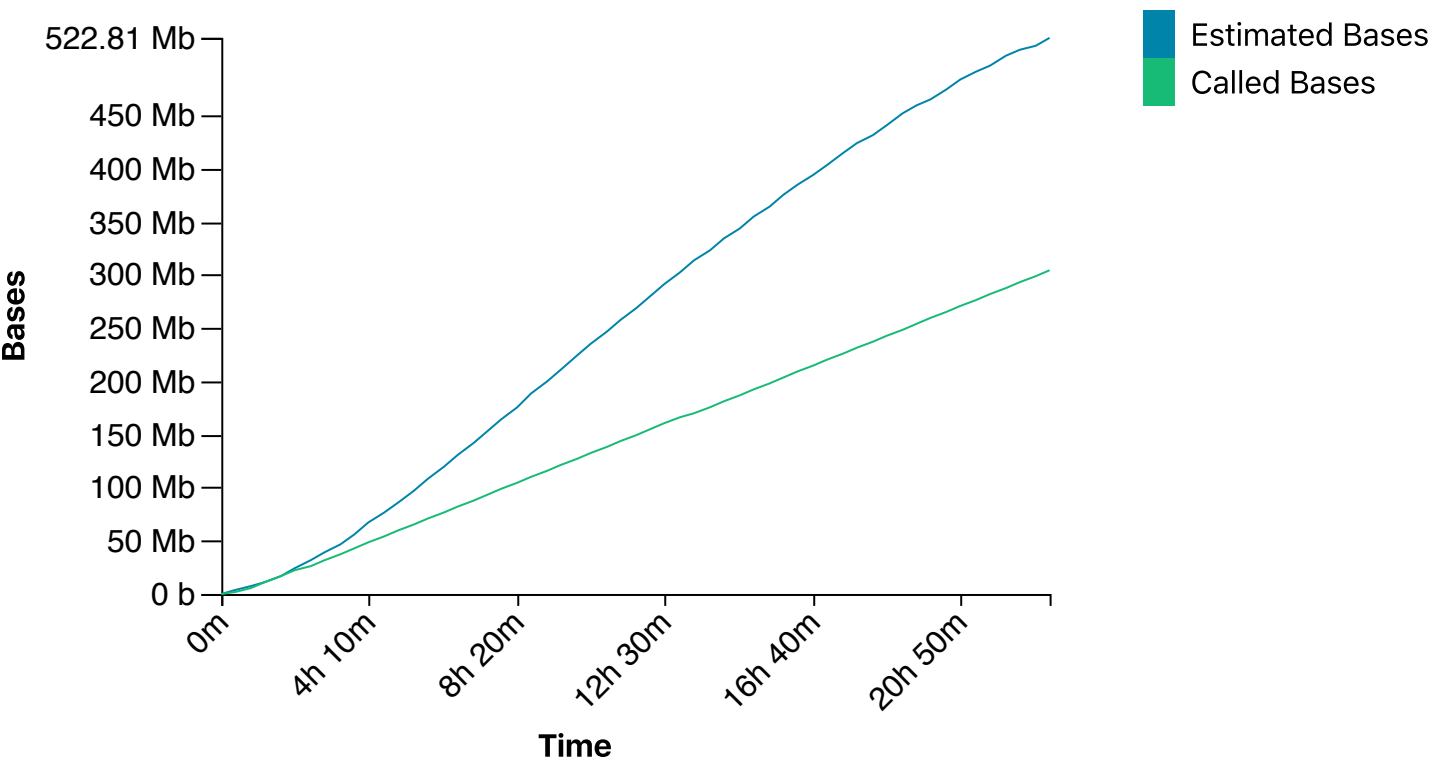


Read Length Histogram Estimated Bases

Estimated N50: 325 b

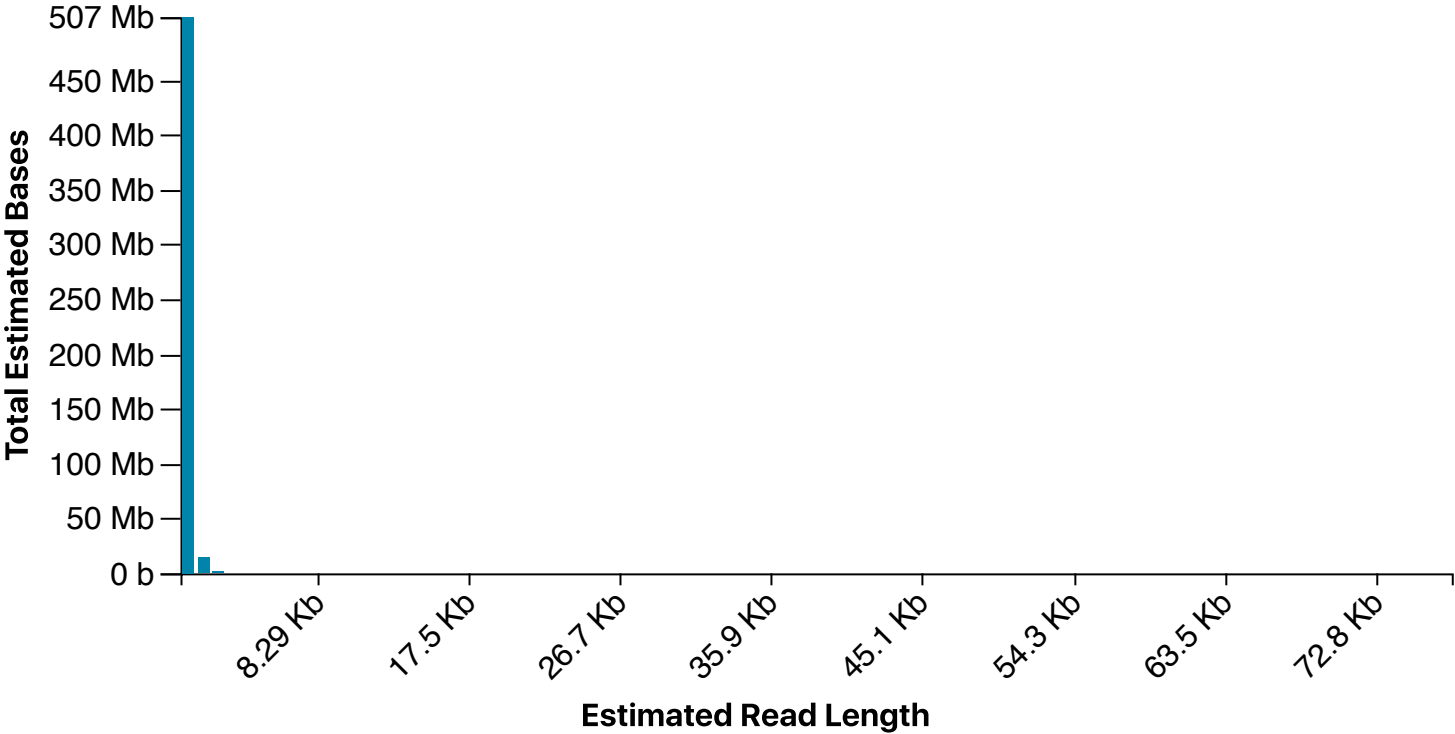


Read Length Histogram Basecalled Bases


Estimated N50: 305 b

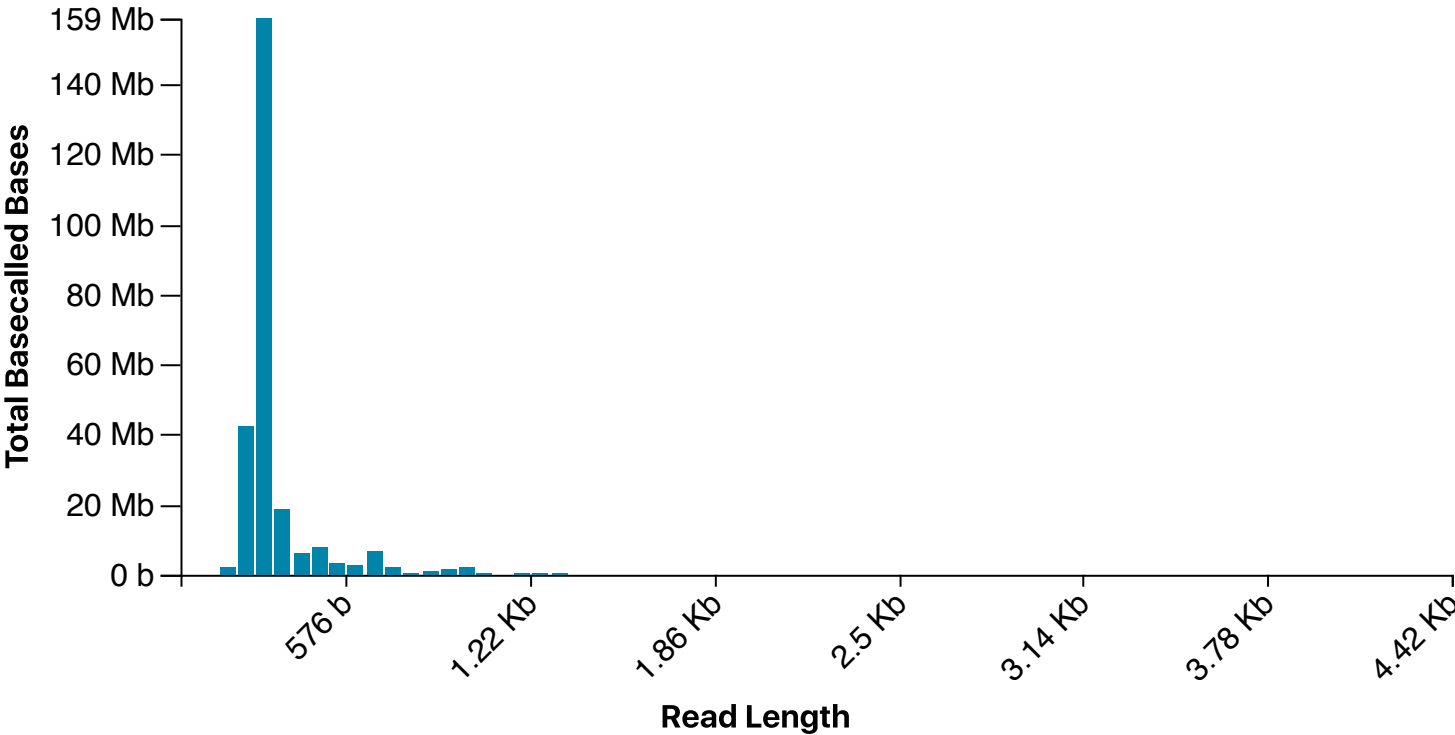


Duty Time Grouped

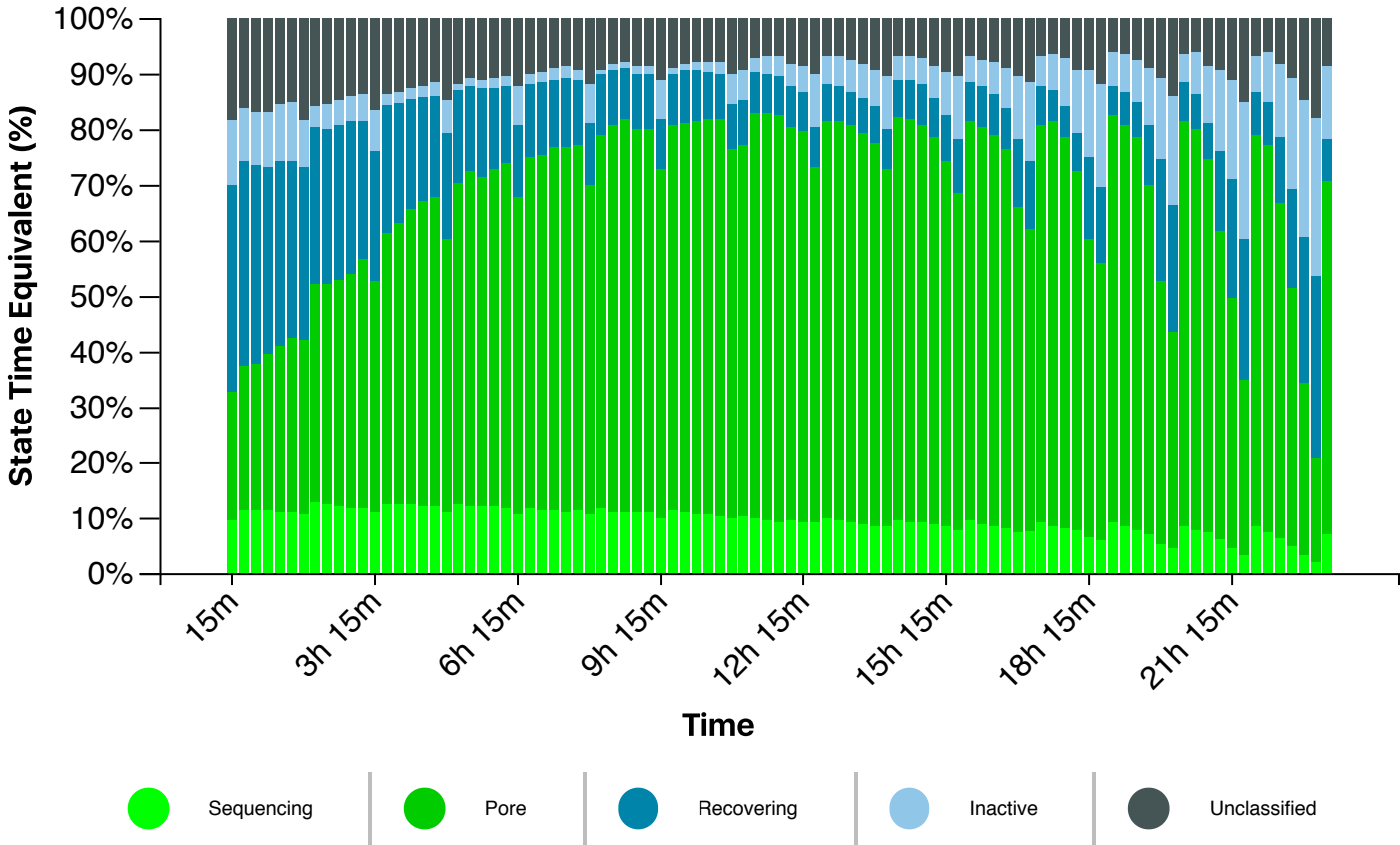


Duty time Categorised

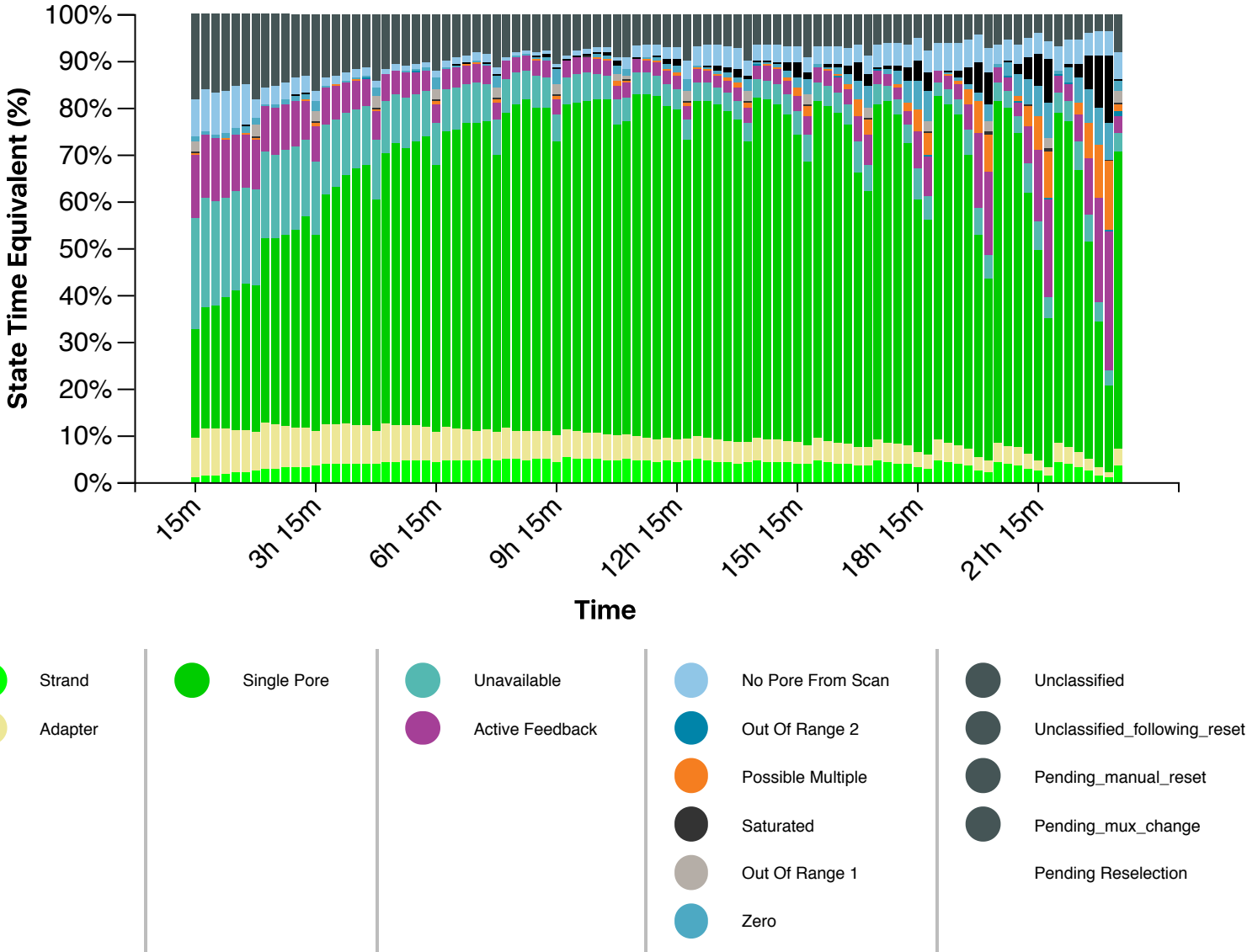


Mux Scan Grouped

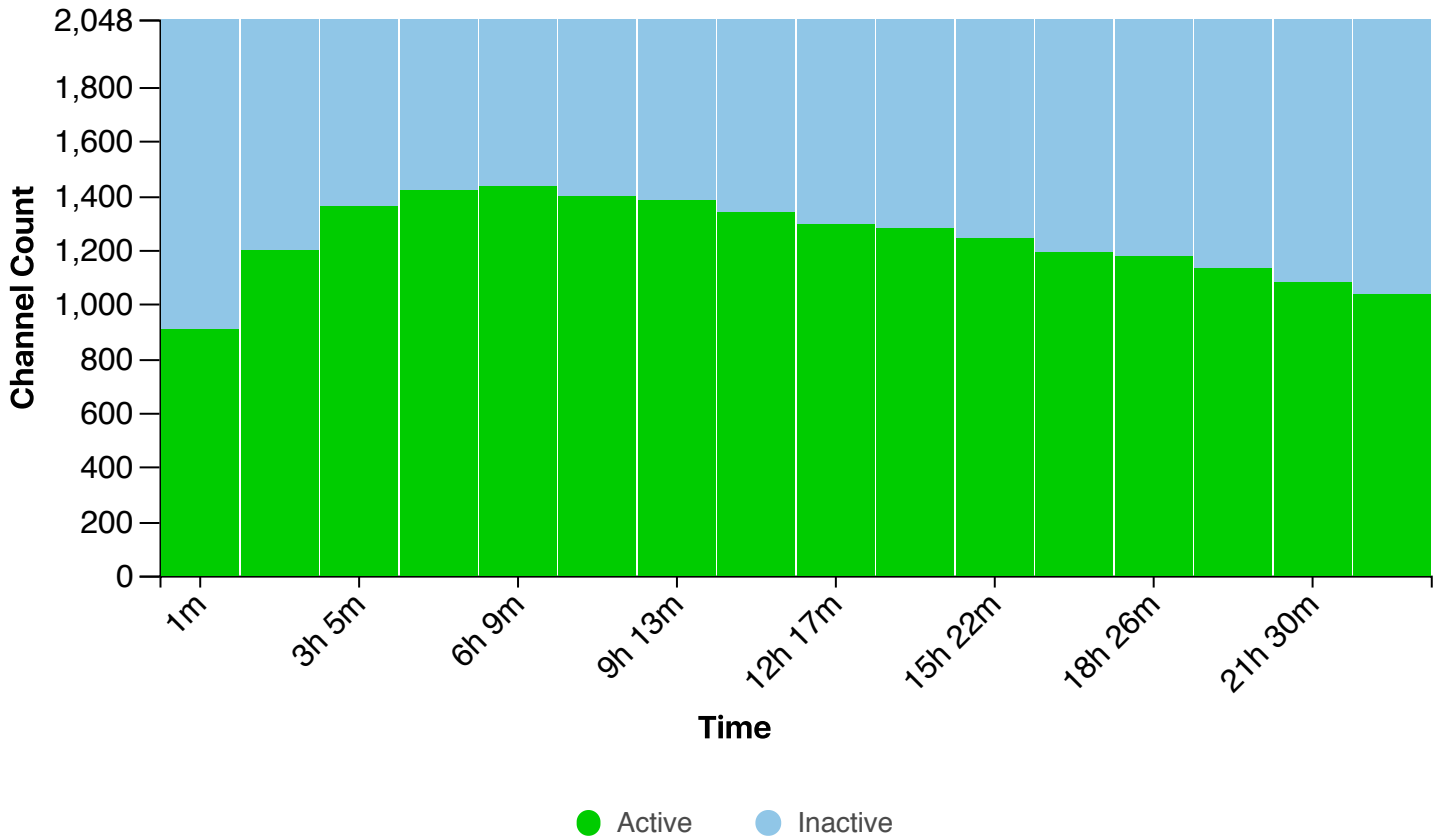


Mux Scan Categorised

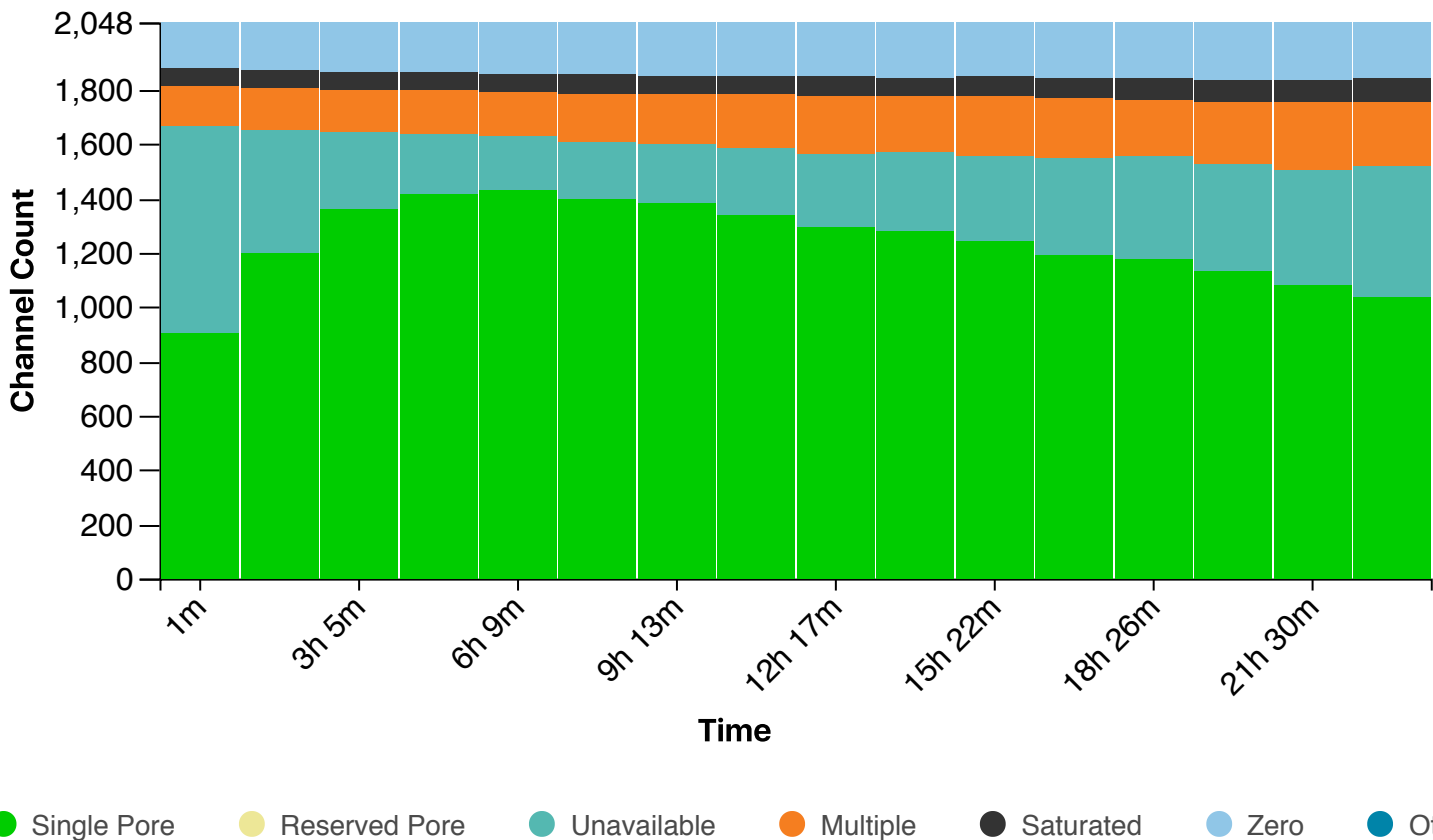


Temperature History

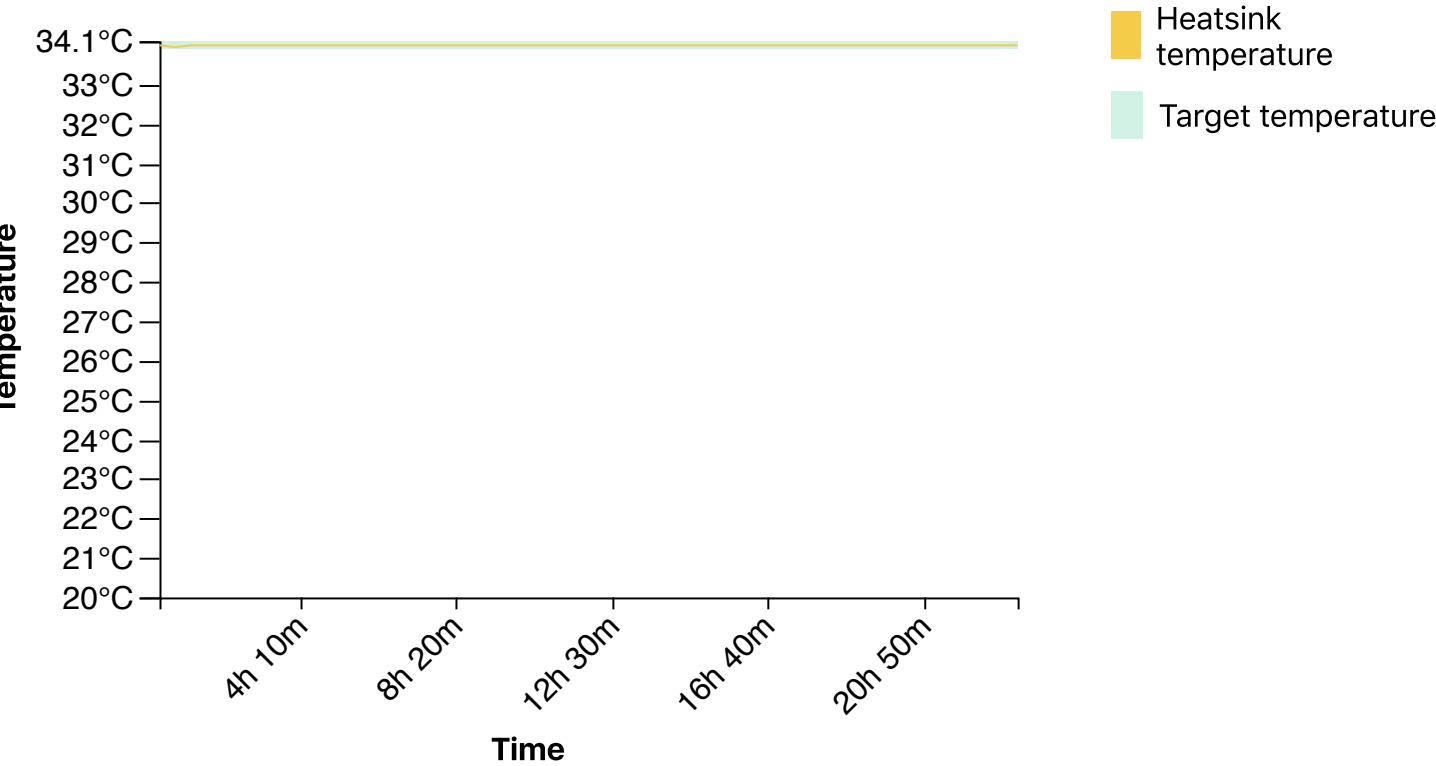


Bias Voltage History

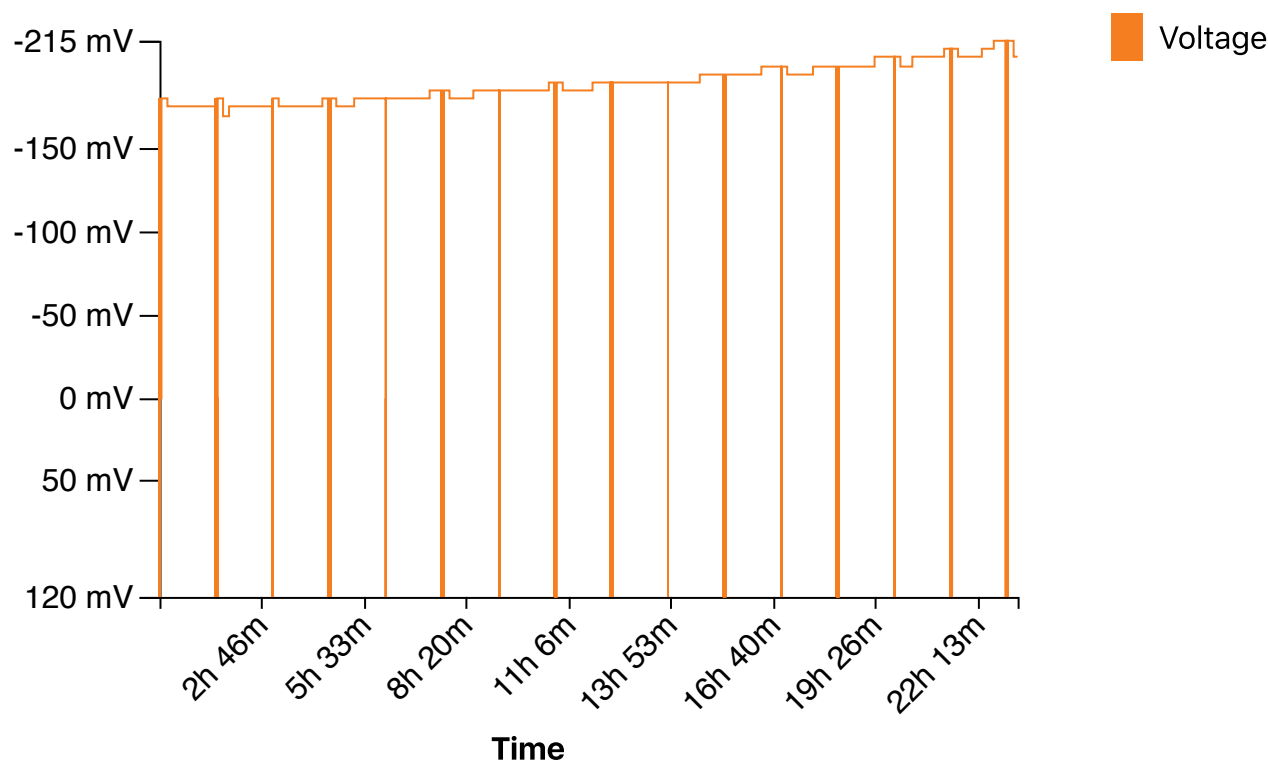


Translocation Speed

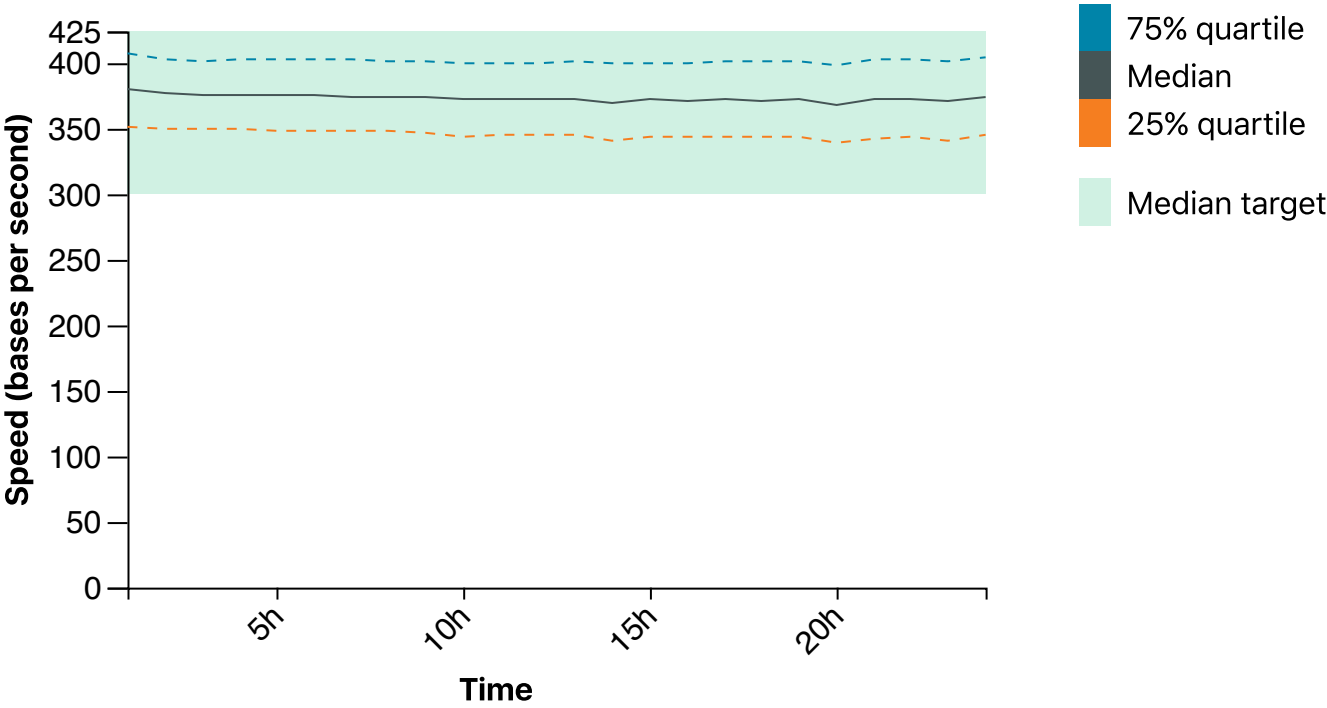

QScore

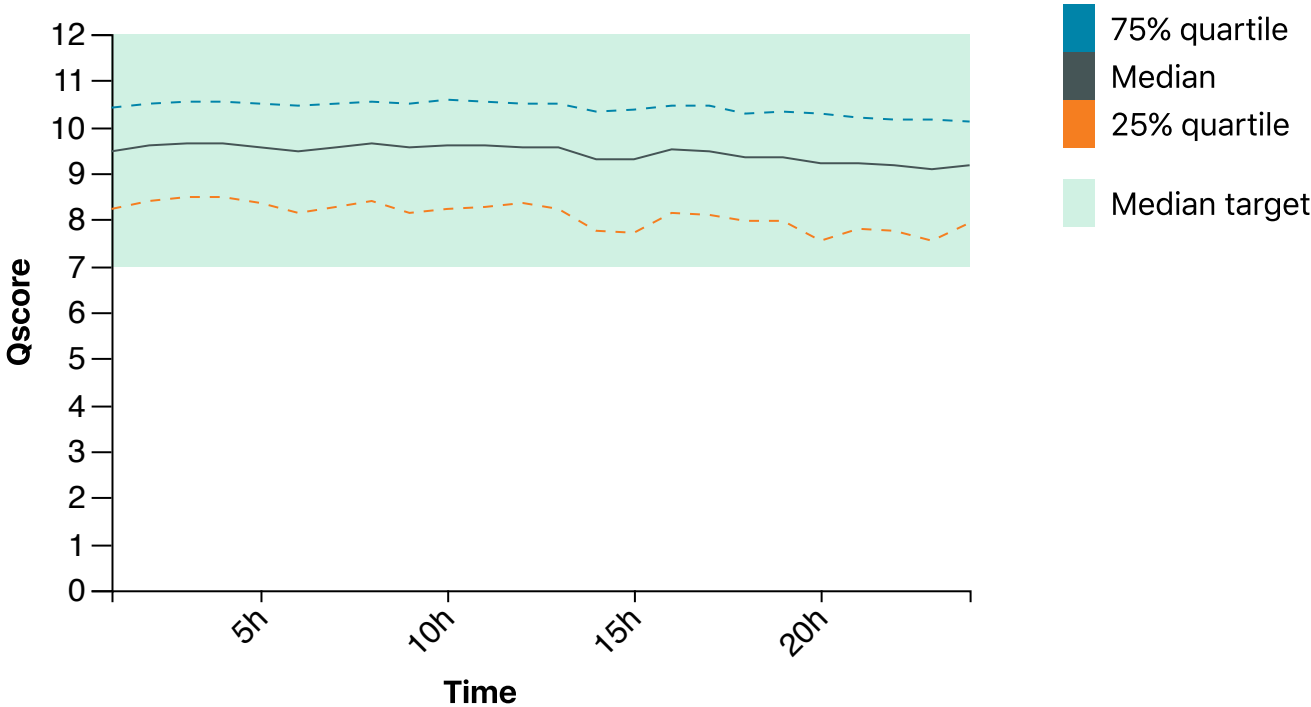

Disk Write Performance

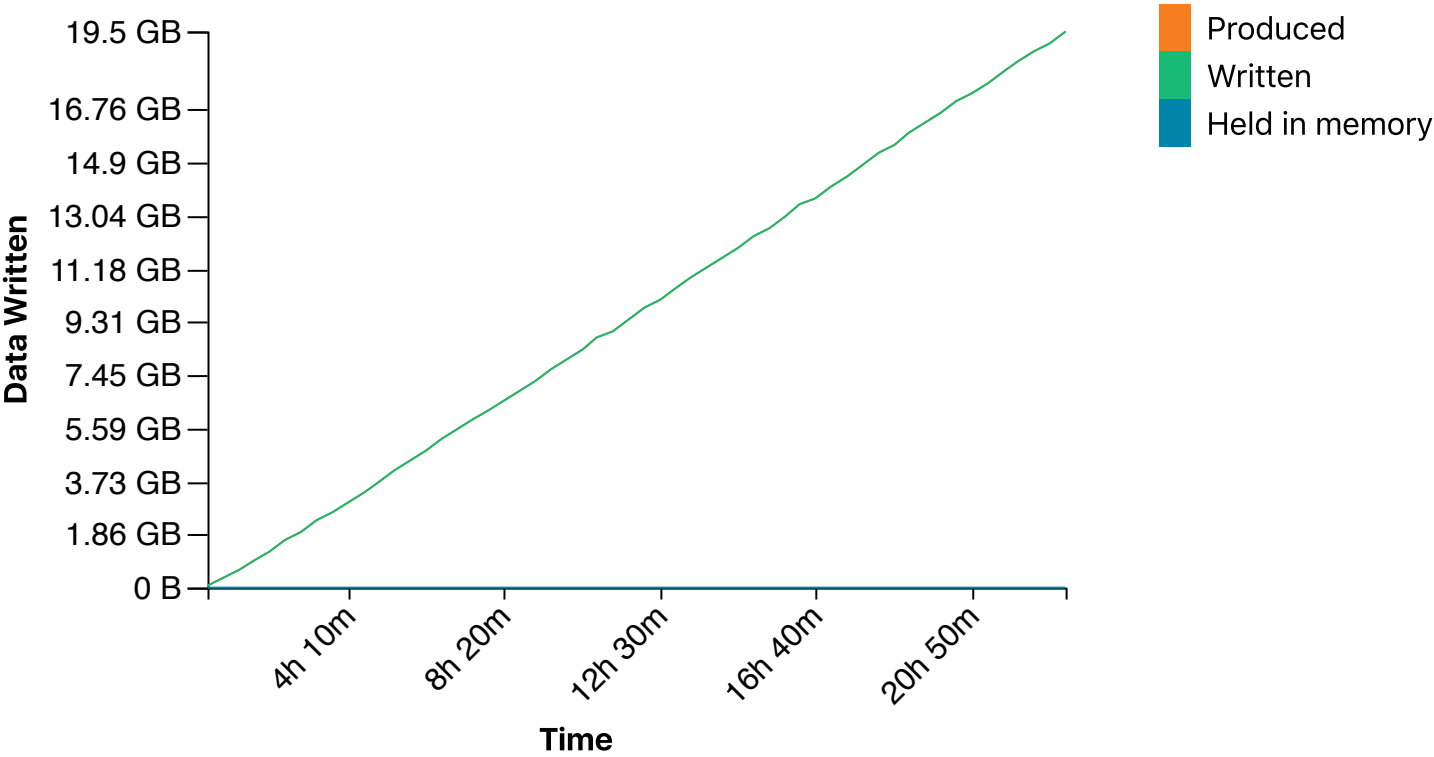

Run Debug Messages

- Flow cell FAN19575 has 1038 pores available for sequencing. Starting sequencing with 479 pores September 22, 13:49
- Performing Mux Scan September 22, 13:48
- Flow cell FAN19575 has 1082 pores available for sequencing. Starting sequencing with 485 pores September 22, 12:17
- Performing Mux Scan September 22, 12:16
- Flow cell FAN19575 has 1133 pores available for sequencing. Starting sequencing with 489 pores September 22, 10:45
- Performing Mux Scan September 22, 10:43
- Flow cell FAN19575 has 1176 pores available for sequencing. Starting sequencing with 481 pores September 22, 09:13
- Performing Mux Scan September 22, 09:11
- Flow cell FAN19575 has 1191 pores available for sequencing. Starting sequencing with 486 pores September 22, 07:41
- Performing Mux Scan September 22, 07:39
- Flow cell FAN19575 has 1245 pores available for sequencing. Starting sequencing with 492 pores September 22, 06:09
- Performing Mux Scan September 22, 06:07
- Flow cell FAN19575 has 1280 pores available for sequencing. Starting sequencing with 493 pores September 22, 04:37
- Performing Mux Scan September 22, 04:35
- Flow cell FAN19575 has 1292 pores available for sequencing. Starting sequencing with 489 pores September 22, 03:05
- Performing Mux Scan September 22, 03:03
- Flow cell FAN19575 has 1335 pores available for sequencing. Starting sequencing with 499 pores September 22, 01:33
- Performing Mux Scan September 22, 01:31
- Flow cell FAN19575 has 1384 pores available for sequencing. Starting sequencing with 507 pores September 22, 00:01
- Performing Mux Scan September 21, 23:59
- Flow cell FAN19575 has 1399 pores available for sequencing. Starting sequencing with 508 pores September 21, 22:29
- Performing Mux Scan September 21, 22:27
- Flow cell FAN19575 has 1430 pores available for sequencing. Starting sequencing with 504 pores September 21, 20:57
- Performing Mux Scan September 21, 20:55
- Flow cell FAN19575 has 1422 pores available for sequencing. Starting sequencing with 506 pores September 21, 19:25
- Performing Mux Scan September 21, 19:23
- Flow cell FAN19575 has 1358 pores available for sequencing. Starting sequencing with 503 pores September 21, 17:53
- Performing Mux Scan September 21, 17:51
- Flow cell FAN19575 has 1199 pores available for sequencing. Starting sequencing with 493 pores September 21, 16:21
- Performing Mux Scan September 21, 16:19
- Flow cell FAN19575 has 909 pores available for sequencing. Starting sequencing with 467 pores September 21, 14:49
- Performing Mux Scan September 21, 14:47

- Starting sequencing procedure September 21, 14:47
- Waiting up to 300 seconds for temperature to stabilise at 34.0°C September 21, 14:46
- Disk / has 1631 GB space remaining September 21, 14:46